

MADE IN GERMANY FOR OVER 100 YEARS.  
DEDICATED TO DIGNIFIED CARE.

VÖLKER

## Healthcare Beds





## The perfect combination of choice, style and functionality.

Völker Healthcare Beds are available in a range of designs and colours, combining the essential functions of a medical device with enhanced styling to give a nursing room a home-like ambience for the improved comfort and wellbeing of residents.

Every person has a right to be treated with respect and dignity. Especially people in need of care. We think, develop and manufacture with this in mind.

### Comfortable. The Völker MiS® Micro-stimulation System.



The patented Völker MiS® Micro-stimulation System is based on neuroscience principles, and detects even minor movements. The lying surface then feeds these tactile stimuli back through the mattress to the resident's body. This continuous feedback increases body perception and promotes movement and blood circulation in the skin, improving the resident's comfort when sleeping or lying down, increasing mobility and reducing the risk of bed sores.





## Dignified care through increased mobility: The Vis-a-Vis bed.

Mobilisation can increase the independence and well-being of residents in long-term care and rehabilitation facilities and encourage them to take part in everyday activities. The Völker Vis-a-Vis bed enables

residents to be brought safely into an upright position. From this position, a single nurse can safely and easily move the resident forwards until their feet are firmly planted on the floor.

## Patented. The Völker siderail system.



The patented Völker siderail system is characterised by its exceptional versatility. It consists of separate elements that are split on both sides and can be adjusted to either one or two levels. This ensures that residents feel adequately protected from accidentally rolling out of the bed. As the siderails are highly versatile and robust, they can also assist mobilisation.

The side elements are directly integrated into the lying surface. The rails at the head and foot move with the bed as it is adjusted.



## Increasing safety by minimising the risk of falls: The 5384 Kepler low-height bed.

The 5384 Kepler provides maximum protection for residents who have an increased risk of falling out of bed. Falls can lead to loss of independence, as well as having other negative implications for both the resident and the care facility.


This versatile low-height bed can be lowered to 23 cm\* and raised to a stable

working height of up to 80 cm\* to ensure a safe and ergonomic working position for caregivers. The bed's close proximity to the floor when at its lowest position can give residents who are prone to falls a greater feeling of safety and may also reduce the occurrence and severity of falls.

\* Dependent on the configuration chosen



→ The ergonomic design of the low-height bed combines light weight, high stability and a safe working load of 230 kg.

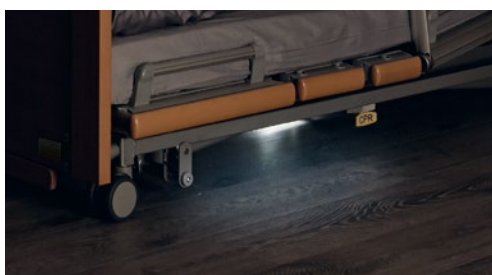

→ The under-bed light improves orientation at night. Controlled using the hand control unit (optional).





## The 5384 Kepler low-height bed at a glance:

- The patented Bed Height Memorisation allows residents and caregivers to easily select a preferred height and position.
- With the EasyChair feature, residents can be brought effortlessly into a seated position.
- Automated regress positioning ensures safer mobility for residents.
- A bed extension and linen holder are optionally available.
- A subtle night light under the bed enhances safety for both residents and nurses (optional).
- Optional double castors increase manoeuvrability and provide automatic braking for enhanced safety when the bed is in its lowest position.

## All Völker Healthcare Beds at a glance.



Model 5384 Kepler




Model 3080



Model 2080



Model Vis-a-Vis

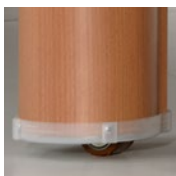

→ On models 2080 and 3080, the castors can be retracted.

## Head and Foot Board designs.

For models 2080, 3080 and 5384 Kepler



Design S



Design FS

For model 3080



Design MA



Design MAH  
(lower foot board/Design MA  
head board)



Design MB

For model Vis-a-Vis



Design S

## Play it safe. With Völker Service.

Völker Service maximises the return on your investment in Völker beds. Our Service team will also take responsibility for a large amount of the work and requirements that arise from the German Medical Devices Act and other legislation.

All of our service technicians are qualified and

experienced testers who are very familiar with the medical devices they inspect. They will test and service your beds meticulously, observing all of the legal requirements.

Any defects are usually fixed immediately and therefore at optimum cost.

Völker Service additionally compiles detailed legal documentation and checks all records to ensure that they are complete.

We would be happy to discuss which Völker Service options would best suit your needs. Simply get in touch.

**VÖLKER**

MADE IN GERMANY FOR OVER 100 YEARS.  
DEDICATED TO DIGNIFIED CARE.

Völker reserves the right to make changes to the design, specifications and models without prior notice. Products may vary slightly from the colours shown here. Subject to printing errors. All information is non-binding.

Not all products/options are available in all countries. For further information on our products or services, please contact your local Völker representative or visit our website.

©2018 Völker | All rights reserved

**VÖLKER**

[www.voelker.de](http://www.voelker.de)



**HEALTHCAREBEDS**  
**AUSTRALIA**

*[info@healthcarebedsaust.com.au](mailto:info@healthcarebedsaust.com.au)*